



**DRG® Benzodiazepine (RAP-4916)**



Revised 15 Nov. 2011 rm (Vers. 4.1)



*This kit is intended for Research Use Only.*

*Not for use in diagnostic procedures.*

*Please use only the valid version of the package insert provided with the kit.*

**INTENDED USE**

The Benzodiazepine (BZO) Test Device is a lateral flow, one-step immunoassay for detection of Benzodiazepine in human urine.

**SUMMARY**

The BZO-Test Device is a competitive immunoassay for detection of benzodiazepine. The detection limit of the test is 300 ng/mL.

The Test Device use highly-specific antigen-antibody reactions for the detection of drugs in human urine.

The test is interpreted visually and provides a qualitative result. In case there is no drug in the urine (non-consumer) two red lines (test- and control-line) appear at the reaction zone of the device. If there is the according drug in the urine sample (consumer) only one control-line appears at the reaction zone.

**1.1 Relevant Information**

This assay provides only a preliminary analytical test result. A more specific alternative chemical method must be used in order to obtain a confirmed analytical result. Gas chromatography/mass spectrometry (GC/MS) has been established as the preferred confirmatory method by the Substance Abuse and Mental Services Administration (SAMSHA). Clinical consideration and professional judgment should be applied to any drug test result, particularly when preliminary positive results are indicated.

**Times of detection of various drugs of abuse in urine (2)**

Drug of abuse/psychotropics	Times of detection
Benzodiazepines	3 days for therapeutic doses, up to 4 to 6 weeks for chronic users


**SET-UP OF THE TEST DEVICE**

The plastic case of the test cassette encloses one test strip. At the right end of the strip there is the sample well and at the left part there is the opening of the reaction zone. At the reaction zone you find the test (T) and the control (C) zone. The labelling next to the reaction zone marks the test and control zone. Because the strip encloses in a plastic case you can only imagine its position by the openings in the case.

## DRG® Benzodiazepine (RAP-4916)



Revised 15 Nov. 2011 rm (Vers. 4.1)

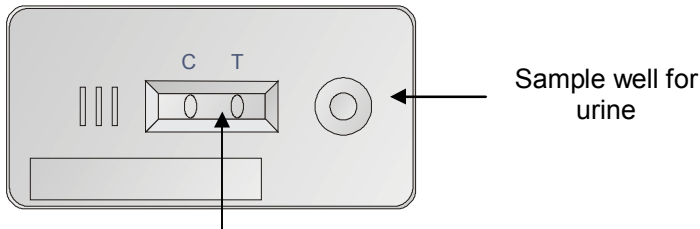

Reaction zone with the test- (T) and the control- (C) zone (with ellipses marked)

**TEST PRINCIPLE**


This test is a competitive immunoassays in which immobilised drug conjugate from the test competes with free drug which may be present in urine for limited antibody binding sites. The test device contains a membrane strip which is pre-coated with immobilized drug conjugate as antigen in the test zone (T). Dye-labelled anti- drug -antibodies are placed in the red coloured conjugate pad below the sample well.

After soaking the membrane with urine by dropping urine in the sample well the antibodies move upwards by capillary action and get to the test zone (T). If there is no drug in the urine the labelled antibodies attach to the immobilized drug conjugate and a visible line is formed. Therefore, a line in the T-zone indicates that no drug is present in the urine (negative result) or that the drug concentration is below the cut-off.

If drug molecules are present in the urine, they compete with the immobilized drug conjugate in the test zone for the limited antibody binding sites. With increasing concentrations of the drug in the sample the binding of the antibody is more and more inhibited and the color of the test line becomes weaker. When the amount of drug is equal or more than the cut-off, it will prevent the binding of the antibody to the drug conjugate of the test zone and the line gradually vanishes. Therefore, the absence of a colored band in the test zone (T) indicates a positive result.

A control line with a different antigen/antibody reaction is also added to the immunochromatographic membrane strip at the control zone (C) to indicate that the test has been performed properly. The presence of this control line serves as verification that sufficient volume has been added, and that proper flow was obtained. The control line should always appear, regardless of the presence of the examined drug.

This means that negative urine will produce two colored lines (non-consumer), where as positive urine will produce only one colored line in the reaction zone (consumer).

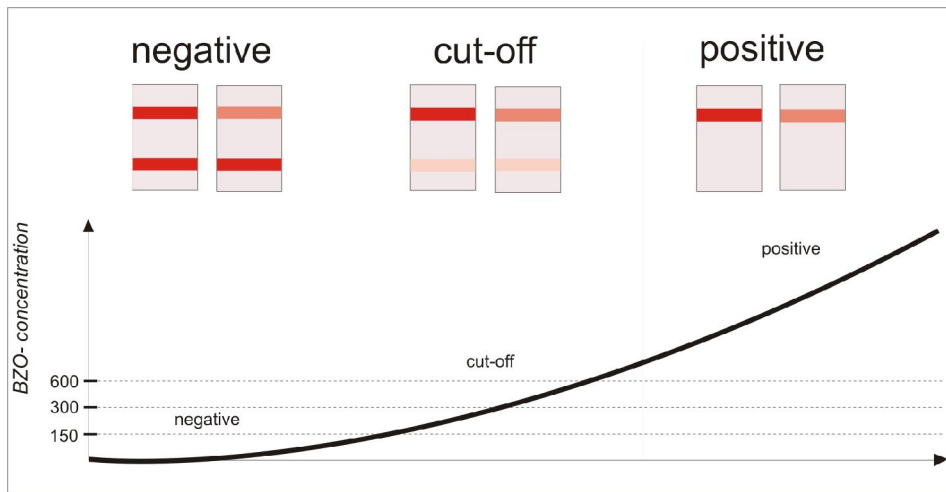

### STORAGE AND STABILITY

The test is to be stored refrigerated or at room temperature (2 °C – 30 °C) in the sealed pouch for the duration of the shelf life.

### PRECAUTIONS

- For professional research use only!
- Use only once!
- Do not use more than the required amount of liquid.
- Avoid cross-contamination of urine samples by using a new specimen collection container and specimen pipette for each urine sample.
- Do not touch the reaction zone of the device to avoid contamination.
- Do not spill the samples into the reaction zone.
- Use only urine as liquid and no other one instead.
- Urine specimens may be potentially infectious. Proper handling and disposal methods should be established.
- Do not use the cassette after expiration date.
- Use test right after unwrapping because it is humidity-sensitive!
- Therefore do not use the test after damage of the packaging foil!
- Please be aware of the developing time of the test before evaluation.
- Please take the specificity and the cross reactivity into account for evaluation.
- Store and transport the test device always at the stated temperature.



## DRG<sup>®</sup> Benzodiazepine (RAP-4916)



Revised 15 Nov. 2011 rm (Vers. 4.1)



- The used test should be discarded according to local regulations.

### REAGENTS AND MATERIALS SUPPLIED

- 50 Individually wrapped test device
- Disposable pipette
- One instruction sheet

### MATERIAL REQUIRED BUT NOT PROVIDED

- Specimen collection container
- Timer

### SPECIMEN COLLECTION AND HANDLING

Collect the urine specimen in a clean and dry collection container. Ensure that a sufficient quantity of the specimen is collected.

This drug test device is formulated for use with urine specimens. Fresh urine does not require any special handling or pretreatment. The test of the urine specimens should be performed as soon as possible after the specimen collection, preferably during the same day. The specimen may be refrigerated at 2 °C - 8 °C for two days, or frozen at -20 °C for a longer period of time. Specimens that have been refrigerated must be equilibrated to room temperature prior to testing. Specimens previously frozen must be thawed, equilibrated to room temperature, and mixed thoroughly prior to testing.

#### Note:

Urine specimens and all materials coming in contact with them should be handled and disposed of as if capable of transmitting infection. Avoid contact with skin by wearing gloves and proper laboratory attire.

### LIMITATIONS OF PROCEDURE

- The assay is designed for use with human urine only. Due to too low specific weight (e.g. absence of ions) in other mediums it might be possible, that false results occur.
- A positive result with any of the tests indicates the presence of a drug/metabolite in urine only, and does not indicate or measure intoxication.
- There is a possibility that technical or procedural errors as well as other substances and factors not listed may interfere with the test and cause false results. See chapter “Specificity” for lists of substances that will produce positive results, or that do not interfere with test performance.
- If it is suspected that the samples have been mislabelled or tampered with, a new specimen should be collected and the test should be repeated.

DRG<sup>®</sup> Benzodiazepine (RAP-4916)


Revised 15 Nov. 2011 rm (Vers. 4.1)

USA: **TEST PROCEDURE**

The test should be performed immediately after opening the protective pouch, because the test is humidity-sensitive. Refrigerated cassettes should be brought to room temperature (15 °C – 30 °C) before opening to avoid condensation of moisture on the test. Refrigerated Patient's samples and controls should be brought to room temperature (15 °C – 30 °C) prior to testing. Ensure that a sufficient quantity of the specimen has been collected.

1. Open the pouch and remove the cassette. Once opened, the cassette must be used immediately.
2. Draw the urine sample to the line marked on the pipette.
3. Dispense 3 drops (120 - 150 µL) into the sample well. Use a separate pipette and device for each sample or control. Start the timer.

Caution: If the urine directly touches the reaction zone, the test gets invalid!

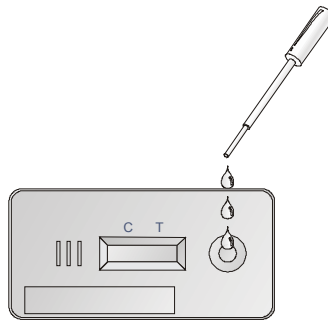


4. At the end of five minutes read the result.

### INTERPRETATION OF RESULTS

For interpretation the reaction zone is regarded. One or two red lines will appear there.


#### Negative Result:

Two colored lines appear. The line in the test zone (T) is the drug line; the line in the control zone (C) is the control line, which is used to indicate proper performance of the strip. The color intensity of the test line may be weaker or stronger than that of the control line.

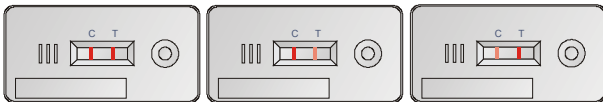

#### Positive Result:

Only one red colored line appears in the control zone (C). The absence of a red test line indicates a positive result. Keep the sample for further confirmation of the result.

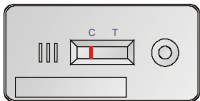


#### Note:

A very faint line in the test zone indicates that the drug in the sample is near the cut-off level of the test. These samples should be re-tested or confirmed with a more specific method before a positive determination is made.

#### Invalid Result:

No line appears in the control zone (C). Under no circumstances should a positive sample be identified until the control line forms. Repeat the assay with a new test device. Paying special attention to the required amount of liquid and reaction time. If problems persist, contact your manufacturer.

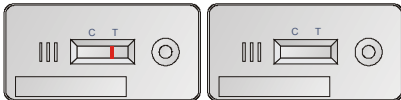

### QUALITY CONTROL

Good laboratory practice recommends the use of control materials to ensure proper kit performance. Quality control specimens are available from commercial sources. For testing the positive and negative controls are used in the same way as urine specimen.

<b>COMPOUND</b>	<b>CONCENTRATION (ng/mL)</b>
Oxazepam	300
a-Hydroxyalprazolam	2,500
Alprazolam	250
Bromazepam	625
Chlordiazepoxide	3,750
Clobazam	100
Clonazepam	2,500
Clorazepate	5,000
Delorazepam	2,500
Desalkflurazepam	500
Diazepam	300
Estazolam	5,000
Flunitrazepam	375
Lorazepam	1,500
Lormetazepam	2,000
Midazolam	10,000
Nitrazepam	37,500
Norchlordiazepoxide	250
Nordiazepam	750
Temazepam	150
Triazolam	1,875

## DRG® Benzodiazepine (RAP-4916)



Revised 15 Nov. 2011 rm (Vers. 4.1)

USA: 

All following listed compounds reacted negative up to a concentration of 100 µg/mL.

Acetamidophene	Guaiacol Glyceryl Ether
Acetone	Hemoglobin
Albumin	Imipramine
Amitriptyline	(+/-)-Isoproterenol
Ampicillin	Lidocaine
Aspartame	(+)-Naproxen
Aspirin	Oxalic Acid
Atropine	Penicillin-G
Benzocaine	Pheniramine
Bilirubin	Phenothiazine
Caffeine	Phenylethylamine
Chloroquine	Procaine
(+/-)-Chlorpheniramine	Quinidine
Chlorpheniramine	Ranitidine
Creatine	Riboflavine
Dexbrompheniramine	Sodium Chloride
Dextromethorphan	Sulindac
4-Dimethylaminoantipyrine	Thioridazine
Dopamine	Trifluoperazine
Erythromycin	Trimethobenzamide
Ethanol	Tyramine
Furosemide	Vitamin C
Glucose	

**LITERATURE**

1. Baselt, R.C. Disposition of Toxic Drugs and Chemicals in Man, Biomedical Publications, 8. Edition, 2008
2. Urine Testing for Drugs of Abuse. National Institute on Drug Abuse (NIDA), Research Monograph 73, 1986
3. Thomas L. eds., Labor und Diagnose, 6. ed., TH-Books Verlagsgesellschaft, Frankfurt, 2005
4. Fed. Register, Department of Health and Human Services, Mandatory Guidelines for Federal Workplace Drug Testing Programs, 53, 69, 11970, 1988
5. McBay, A.J. Clin. Chem. 33, 33B-40B, 1987.
6. Hofmann F.E.: A Handbook on Drug and Alcohol Abuse: The Biomedical Aspects, New York. Oxford University Press, 1983.

Version 2011-11-08~rm